

Instructions for use


Product Name: Orthodontics Elastics
Type: Ligature Elastics, Separating Elastics, Intro-oral Elastics
The product used for: For orthodontic treatment only.
Where to use: This product is used in orthodontic/dental clinic.

Clinic Practice Advice:


1. Choose the exact right products according to your need.
2. Before using, the product is suggested to be cleaned by 75% alcohol, Soaked in alcohol in 20 minutes. There should be no change on the color and the shape of the products after cleaned and autoclaved, as well as no impact to the treatment.
3. The tooth surface should be cleaned , dried and isolated from wet circumstance as well. Following the advice of professional textbooks and training course, the orthodontic elastics use together with the orthodontic brackets, wires, etc.

Key Points to Success:

1. For the Separating Elastic, suggest to use the Separator placing pliers to make it easier to put the modules between the tooth. Keep 48hours and then take off. The purpose of this product is to make enough space to ensure the bands can be put onto the teeth.

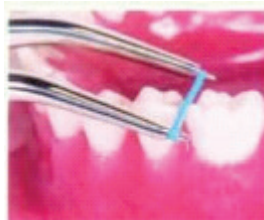

2. For the Ligature Elastic, the Needle Holder plier is recommended to use to fix the wires to the brackets. the usage time should be within one week.

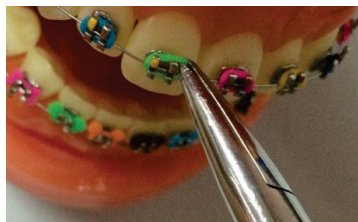

Warning:

1. The product contain Rubber, rubber allergies case should be prohibited to use!
2. Single use only! Do not re-use! Repeated use can cause mechanical performance failure and cause cross-infection.
3. For professional use only. Mixed up or wrongly placement of this products might lead to unexpected problem in clinic .

Storage

Be stored in dry and room temperature. Avoid the contact of corrosive airs.